

Handheld Pulse Oximeter

Model: SA310



- Big 2.4" color TFT display
- Personal SpO2, Pulse Rate, Memory Interval thresholds setting
- Shielded Design blocks Ambient Light
- Biocompatibility & Anti-Allergic Design
- Alarm function
Device delivers audible warning when measurement result exceeds the SpO2/Pulse Rate thresholds previously set.
- Silent mode
Disables the beep sound during the measurement by pressing the "Sound" button.
- Up to 999 Memories
Could be programmed with 3 intervals: 10 / 30 / 60 seconds
- Bluetooth® Transmission
- Pulse Strength Indicator
- Dual power source
- Protection pouch



%SpO2 & Pulse Rate




2.4" Color Display



High/Low SpO2, PR & Memory Interval Setting



Shielded Design blocks Ambient Light

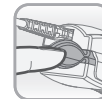

Biocompatibility & Anti-Allergic Design



Alarm & Sound Off



Up to 999 memories



Pulse Strength Indication



Battery and AC Adaptor



Protection Pouch



Bluetooth® Transmission

The Bluetooth® word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by Rossmax International Ltd. is under license. Other trademarks and trade names are those of their respective owners.

Optional Accessory:



Adult Probe PA100



Neonatal Probe PB100



Pediatric Probe PC100



Disposable Probe PD100



Tube Probe PF100

Pulse Oximeter

- Ideal for oxygen therapy assessment
- Ideal for sleep apnea screening
- Ideal for health management
- Ideal for spot-check & continuous monitoring

Specification

Measurement Range | SpO2: 35%~100%; PR: 30~250 bpm

Precision | SpO2: 70%~100%; ±2%; 35%~69% (unspecified); Pulse rate: 30~250bpm; ±3

SpO2 Accessory | Adult Probe Included, Pediatric / Neonatal Probe Optional

Power Supply | 4 x "AA" Alkaline / AC Adaptor

Model	Qty per carton	Carton volume
SA310	10 pcs	0.055 cbm/ctn

- Storage temperature: -25°C - 70°C (-13°F - 158°F)
- Relative humidity: 15% - 90% (no condensing)
- Atmospheric Pressure: 700 to 1060 hPa